GEBRAUCHSANWEISUNG
INSTRUCTION LEAFLET
INSTRUCTION D'UTILISATION
INSTRUCCIONES DE EMPLEO
ISTRUZIONI PER L'USO
GEBRUIKSAANWIJZING
РУКОВОДСТВО К ПРИМЕНЕНИЮ

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MANY THANKS FOR YOUR CONFIDENCE ! A OPERATION OF THE INSTUMENT

We are happy that you have decided to purchase a **Promed** product. You are now the owner of a high-quality product for body care and cosmetic applications. Promed is one of the leading companies with fully technically developed devices for skin and nail care, proven and long-standing.

Devices which make it easy to treat oneself or can be used for cosmetic care problems.

Expert advice is available everywhere where **Promed** articles are sold, or you can simply ask us directly. We can name you a competent consultant.

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Lindenweg 11, D-82490 Farchant, Tel: 08821-9621-0,

Fax: 08821-9621-21, eMail: info@promed.de

Further information can be found on our homepage

www.promed.de

- The device will only work connected to the mains. Therefore, before starting the device connect the mains adapter to the socket.
- For devices with a mains cable, connect the adaptor to the mains socket before operating to use with mains power. For battery-operated devices, you should charge up the battery before operating, using the wall plug transformer according to the instructions.
- The device is now ready.
- Switch the device on by turning the speed control and set it to the required speed using the same control.

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

B TIPS ON CARE

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

Skin and nail care is not merely a matter that is indicated when suffering correspondingly, it is also a cosmetic matter.

Start with cutting and forming toenails. They absolutely must be filed straight in order to hinder nails growing into the skin. Micro-fine filing using the sapphire disc protects against ingrown nails.

Troublesome calluses always develop in the areas where the foot is over-stressed, mostly the heel, the front balls and the side of the big toe. Treat the callus using the sapphire cone or even better with the drainage grinder. Repeat this procedure regularly. Under no circumstances should you take radical measures when removing the callus.

Ingrown toenails can – as needed – be loosened gently using the small cutter. Thick lignified nails can be ground off slowly from the top and filed thinner from the sides.

This takes place without enduring any pain whatsoever and without any danger of injuries. Crumbly nail areas can also be cleaned

using the fine sapphire cutter.

Corns are the painful result of a slow callus process. They are calluses with an externally hardened part which is soft inside. Remove the upper callus with the sapphire cone and then work further using the small sapphire cutter. This already removes the pressure. Deeply sunken corns should not be removed by oneself; they should be removed by a doctor.

If skin and nail care are carried out systematically using the **Pro-med** devices and accessories, regular care becomes a positive experience.

The special abrasive wheels are described more closely on the following pages.

The abrasive wheels are available in various forms and degrees of hardness and can be re-ordered at any time. Ask your specialist dealer.

C APPLICATION OF THE ABRASIVE WHEELS

Abrasive wheel

Shortens and forms the nails right up to the corners. Result: perfectly cared-for nails which no longer tear, break or split.

Abrasive cone

For calluses and corns on feet and hands, also suitable for grinding highly lignified nails. It helps to bring bent nails back to form and eliminates painful pressure.

Callus drainage grinder

Especially suitable for strong, thick calluses on heels and the balls of the feet where there are often several calluses.

Felt polishing cone

The polishing cone to finish all nail care. It gives toenails and fingernails a natural shiny finish. Perfectly smooth nails and nail edges give the care an ideal finish.

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

C APPLICATION OF THE ABRASIVE WHEELS

DE

US

FR

IT

ES


NL

RU

PL


SE

FI

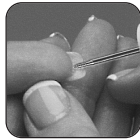Cutter, pointed

For the elimination of difficult corners and edges which thicker grinders cannot reach. (Be careful when using the pointed cutter!)

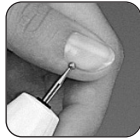
Cutter, round

This is used to grind out ingrown nails. Use it from the sides and cut the ingrown nails out slowly. The round sapphire cutter should not be used to remove nail skin at the moon of the nail as the danger exists that grooves could be cut into the nail moon.

Cutter, long


Grinds ingrown nails free slowly. For instance, it is also possible to remove nail fungus spores easily from under the front of the nail.

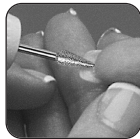Cutter, small cone / stone grinder, cylinder

For the elimination of difficult corners and edges which are not accessible with thicker grinders. (Be careful when using the pointed cutter!)

C APPLICATION OF THE ABRASIVE WHEELS

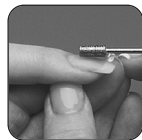
Cutter, barrel / stone cutter, barrel

For grinding thick, lignified nails. Lignified nails can be brought into good form, look better and the pressure eliminated.

Cutter, football

For grinding fungicidal nails; nail fungus spores can be eliminated easily. After removal of the spores, a chemist's tincture should be used to increase the effect. It is also appropriate for grinding hard corn edges to eliminate pressure. Only grind the pressure areas!

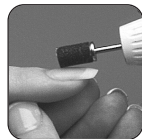Grinding mandrel, plastic + abrasive band, fine/medium/rough

For the preparation of natural and artificial nails, together with gel, acrylate for instance.

Nail sealant

For sealing nail borders and smoothing upper nail surfaces.

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

C APPLICATION OF THE ABRASIVE WHEELS

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

Silicone polisher

Three different grains to polish artificial nails (gel, acrylate etc.).

Polishing brush

To polish up nail setting lotion or ointments against nail fungus.

Soft Polisher

With different layers and materials to polish up artificial nails and natural nails.

D MAINTENANCE AND CLEANING

1. General safety

So as not to impair the safety, reliability and performance of the instrument, expansions, resetting, modifications or repairs should only be carried out by persons authorised to do so by **Promed**.

2. Cleaning the instrument

Only clean the instrument with a soft, damp cloth and a little washing up liquid. Our after-sales service department will be pleased to assist you if you require any further information.

3. Cleaning the sapphire heads

The heads can be cleaned of dust and skin residue with a brush (art. no. 198910).

We also recommend that the heads are disinfected with standard disinfectant spray or 90% alcohol (available from your pharmacist) to prevent the transfer of infections.

E CUSTOMER SERVICE / TECHNICAL DETAILS

Any repair to the monitor may only be carried out by the manufacturer/supplier or by any authorized service center.

The technical details can be found in the description of the specific device

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

F SAFETY HANDLING

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

1. General

The device is radio-shielded according to the EC directive 87/308/EEC. It is designed for network operation at 230 V/50 Hz and is connected to the mains with a transformer.

2. Danger

To reduce the risk of electrocution:

- Do not use the instrument while bathing or washing
- Never place or store the charger and instrument where it can fall or be pulled into a tub, sink or toilet.
- Never reach for an appliance that has fallen into water. Unplug immediately.
- Never submerge the charger in water or any other liquid.


3. Warning

To reduce the risk of burns, electrocution, fire or injury:

- Do not operate the instrument if it has a damaged cord or plug, if it is not working properly or if it is immersed in water.
- Do not excessively bend or pull the cord. Unplug the cord by

grasping the plug - not the cord. Keep the cord away from heated surfaces.

- This product has been designed for manicure-pedicure **ONLY**. Do **NOT** use the instrument for any other purpose.
- Children should be supervised to ensure that they are not playing with the unit.
- This unit is not designed to be used by people (including children) who have limited physical, sensory, or intellectual abilities or who do not have sufficient experience and/or knowledge. This does not apply if these people are being supervised by a person responsible for their safety or receive instructions from them on how to use this unit.

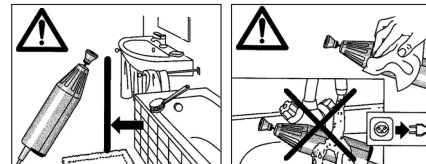G WARRANTY INFORMATION

This instrument has been manufactured carefully and checked thoroughly before leaving the plant. In presentation of warranty card issued for this instrument we therefore give a warranty of 24 month starting on the date of purchase according to the following conditions:

- If in spite of proper use detectable defects concerning material and manufacture should occur and be detected within the period of warranty we will replace all faulty parts of the instrument free of charge, as well as covering any labour charges incurred by the repair.
- The warranty expires if the instrument has not been handled properly (e.g. fracture), after intervention (e.g. opening of the instrument's case) as well as if unauthorized spare parts are used. Parts subject to wear are excluded from the warranty. The devices' engines and moving parts are not part of this guarantee.
- The warranty period starts on the date of purchase. The use of a warranty service does not affect the duration of the warranty. Warranty service must be claimed before the expiry date of the warranty. Complaints which occur after the warranty has expired cannot be taken into account.
- The warranty within the frame of these conditions comes into force only if the date of purchase on the warranty card has been endorsed by the stamp/signature of the dealer.
- In the case of complaint or repair please send the entire instrument with completed warranty card to your appropriate customer service.

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Promed GmbH