<p>(D) Massagesessel Gebrauchsanleitung</p>	2-11	<p>(I) Poltrona per massaggi Istruzioni per l'uso</p>	42-51
<p>(GB) Massage chair Instructions for use</p>	12-21	<p>(TR) Masaj koltuğu Kullanım kılavuzu</p>	52-61
<p>(F) Fauteuil de massage Mode d'emploi</p>	22-31	<p>(RUS) Массажное кресло Инструкция по применению</p>	62-71
<p>(E) Sillón de masaje Instrucciones de uso</p>	32-41	<p>(PL) Fotel masujący Instrukcja obsługi</p>	72-81

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Dear customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of beauty, heat, gentle therapy, blood pressure/diagnosis, weight, massage and air.

1. Intended use.....	12
2. Safety notes.....	12
3. Getting to know your device	14
4. Unit description	15
5. Initial use.....	17
6. Operation.....	18
7. Maintenance and storage	20
8. What if there are problems?.....	21
9. Technical Data.....	21
10. Disposal.....	21
11. Warranty and service	21

Please read these instructions for use carefully and observe the information they contain.

Kind regards,
Your Beurer team

Signs and symbols

The following symbols are used in the instructions for use and/or on the type plate.

Warning!

Warning instruction indicating a risk of injury or damage to health

Caution!

Safety note regarding potential for damage to the device/accessories

Important!

Note on important information

Only use indoors

Do not insert needles!

Cross-reference to another section.

Included in delivery

- Massage chair
- Back cushion
- Power supply unit

1. Intended use

The massage chair is intended solely for massaging the back and legs.

Please observe the safety notes in these instructions for use.

Use the device only according to the methods specified.

The device is suitable only for **private use**, and not for medical or commercial purposes.

2. Safety notes

- Read these instructions for use carefully. Non-observance of the following information may result in personal injury or material damage.
- Store these instructions for use in a safe place and make them accessible to all users.
- Make sure you include these instructions for use when handing over the device to third parties.

Electric shock

Warning!

- As with every electrical device, use this massage chair carefully and cautiously in order to avoid hazards resulting from electric shock.
- Use the device only with the power supply unit provided and at the mains voltage specified on the power supply unit.
- Never use the device if it or any of the accessories show visible damage.
- Do not use the device if it has been accidentally splashed with water.
- Do not use the device during a storm.
- In the event of defects or operational faults, switch off the device immediately. Pull the mains plug out of the socket.
- Never pull on the mains cable to remove the mains plug from the socket.
- Do not clamp or pinch the mains cable.

- Keep the mains cable away from hot surfaces.
- Ensure that the device, switch, power supply unit and cable do not come into contact with water, steam or other liquids.
- Use/operate the device only with dry hands.
- Do not subject the device to any impacts or drop the device.

Repairs

Warning!

- Do not open the device.
- Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. For repairs, please contact Customer Services or an authorised retailer. Failure to do so will invalidate the warranty.
- If the mains plug, mains cable or power supply unit of the device are damaged, they must be replaced by the manufacturer or at an authorised Customer Service point.
- The zip fastening of the massage chair must not be opened. It is only used for technical production reasons.

Risk of fire

Warning!

There may be a risk of fire if the device is used improperly or these instructions for use are ignored!

- Do not use the device together with other electrical devices (e.g. electric blanket).
- Never use the device near petrol or other highly flammable substances.
- Never use the device underneath a cover (such as a blanket, pillow, etc.)

Operation

Warning!

- Do not insert needles or other sharp objects into the chair.
- Do not force or jam any part of your body, in particular your fingers, between the rotating massage heads or their holders inside the device.

This device is exclusively intended to massage the back and legs of humans. It cannot replace medical treatment. Do not use the massager if one or several of the following warnings applies to you.

Do not use the massager

- If you suffer from a medical abnormality or injury to your back (e.g. slipped disc)
- On children
- During pregnancy
- While sleeping
- On animals
- Whilst undertaking any activity where an unexpected reaction could be dangerous
- After the consumption of substances that lead to limited cognitive ability (e.g. pain medication, alcohol)
- Never massage areas that exhibit swellings, infections or eczema.
- Do not use the device in the event of inexplicable stomach ache.
- The device features a switch-on heat function in the massage heads. Persons sensitive to heat must be careful when using this function.

If you are unsure whether the massager is suitable for you, please consult your doctor.

You should consult your doctor before using the massager

- If you suffer from a severe illness or have recently had surgery on your upper body
 - If you have a pacemaker, implant or other aids
 - If you suffer from thrombosis
 - If you suffer from diabetes
 - If you are suffering pain, the cause of which is unknown.
- The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

- This device is not intended for use by people (including children) with restricted physical, sensory or mental skills or a lack of experience and/or a lack of knowledge, unless they are supervised by a person who is responsible for their safety or are instructed by such a person in how to use the device.
- Do not insert any objects into the device apertures or the rotating parts. Ensure that the moving parts can move freely at all times.
- Children should be supervised around the device to ensure they do not play with it.
- Always switch off the device at the control element after use.
- Ensure that the mains cable does not pose a trip hazard.
- Never stand on the device or place any objects on the device.
- Only one person may use the massage chair at any given time.
- Never sit on the arm rest.
- Never load the device with more than 135 kg.

 **Caution!**


- Do not expose the device to direct sunlight.
- Do not use the device at temperatures above 40°C.
- Check the function of the device if it has not been used for a long period of time.

3. Getting to know your device

You can enjoy your own massage at home with the Beurer MC3000 massage chair. The massage chair provides a soothing and relaxing massage.

4. Unit description

Device parts

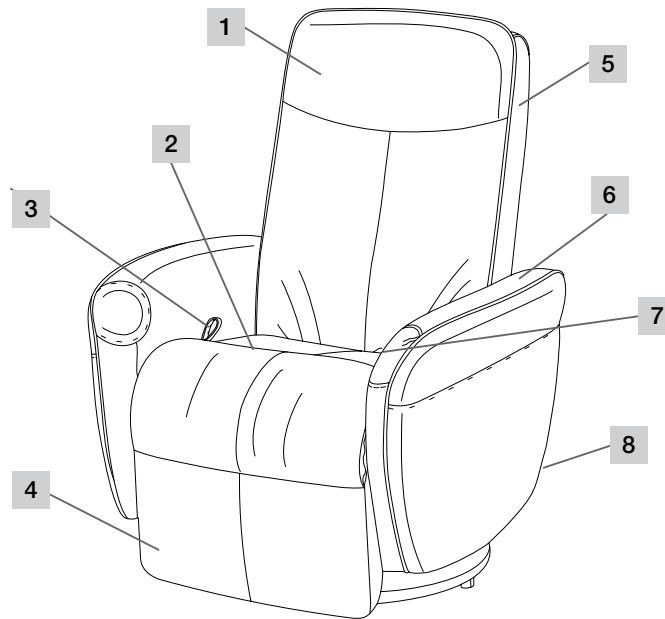1 Back cushion

2 Seat

3 Back rest release

4 Leg rest


5 Back rest

6 Arm rest

7 Control element (interior side of arm rest)

8 Power supply unit with connector plug and mains cable

Control element

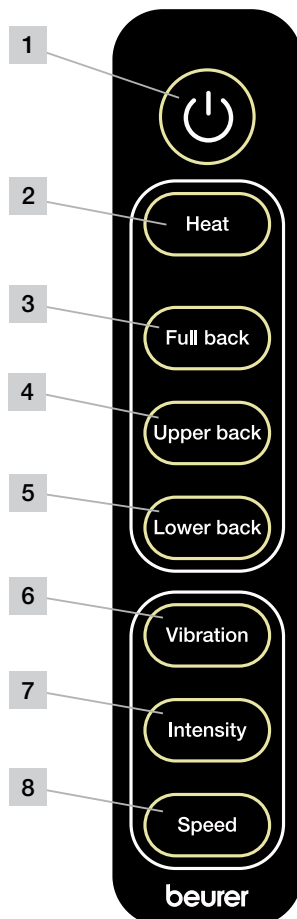1 On/off button

2 Heat, light and heat function

3 Full back

4 Upper back

5 Lower back

6 Vibration

7 Intensity

8 Speed

Control element button functions

On/off button

- > Press the on/off button .
The massage chair is ready for operation.

Massage areas

- The massage can be stopped at any time. To do so, press the on/off button .
The massage mechanism then returns to the start position and switches off.
The currently selected massage area is indicated by the associated LED buttons.


Full back

- > Press the Full Back button .
The 3D Shiatsu massage starts up for the entire back. Total duration: 15 minutes.
The massage mechanism then returns to the start position and switches off.


Upper back

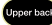- > Press the Upper Back button .
The 3D Shiatsu massage starts up for the upper back. Total duration: 15 minutes.
The massage mechanism then returns to the start position and switches off.

Lower back

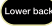- > Press the Lower Back button .
The 3D Shiatsu massage starts up for the lower back. Total duration: 15 minutes.
The massage mechanism then returns to the start position and switches off.

Heat function

The massage heads can be heated.

The heat function can only be switched on when the massage programme is running.

Heat


- > Press the Heat button .
The massage heads are heated.
The LED button lights up.
- > Press the Heat button  again.
The heat function is switched off. The LED button light extinguishes.
When the current massage programme is complete, the heat function automatically switches off.

Settings

Vibration - in-seat and calf-area vibration massage

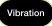- The in-seat and calf-area vibration massage can be switched on.
Three vibration programmes are available:
- > Press the Vibration button .
The vibration massage starts up.
Each time the button is pressed, the next setting is selected.
Pressing the button for a fourth time deactivates the vibration massage.


Intensity - intensity of the massage programme

- There are 3 intensity settings.
- > Press the Intensity button  while the massage programme is active.
The intensity of the massage switches to the next setting.
Each time the button is pressed, the next setting is selected.

Speed

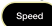- The massage has 3 different speeds.
- > Press the Speed button  while the massage programme is active.
The massage switches to the next speed.
Each time the button is pressed, the next speed is selected.

5. Initial use

Warning!

- Keep packaging material away from children!
Choking hazard.
- Before using the device for the first time, read the 'Safety notes' section.
▷ Safety notes, pages 12 and 13

Caution!

- The massage chair could damage your flooring.
Place a suitable, non-slip underlay under the massage chair.

Installation location

- Install the massage chair on a sturdy and even surface.
- Before using the device for the first time, remove all obstacles in front of and behind the massage chair so that the back rest can be lowered and the leg rest raised freely.


Unpacking and assembly

The massage chair MC3000 is delivered in a cardboard box.

Caution!


Do not use sharp or pointed tools (e.g. a Stanley knife) to remove the packaging.

Step 1:

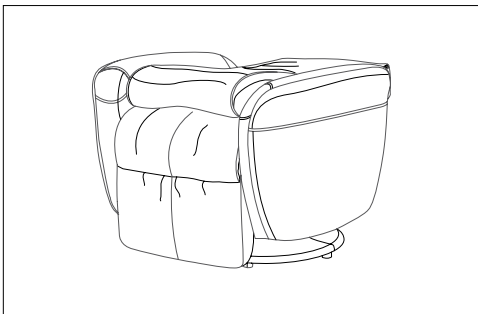
- > Open the cardboard box and hold the massage chair at the designated holding points
- > Remove the chair from the cardboard box
Two people are required for this.
- > Carefully remove all the packaging material

Step 2:

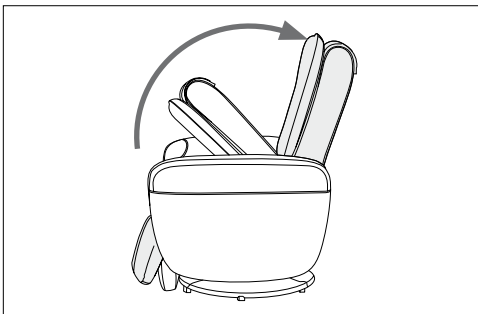
- > Unfold the back rest to the upright position until the locking mechanism audibly engages
The back rest now stands at approx. 115°.

Step 3:

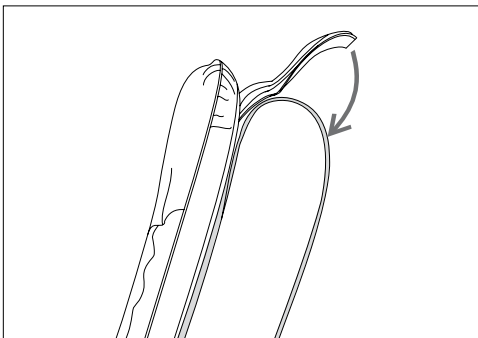- > Lay the back cushion on the massage chair over the back rest
- > Connect the back cushion and the back rest using the zip fastening

Step 4:

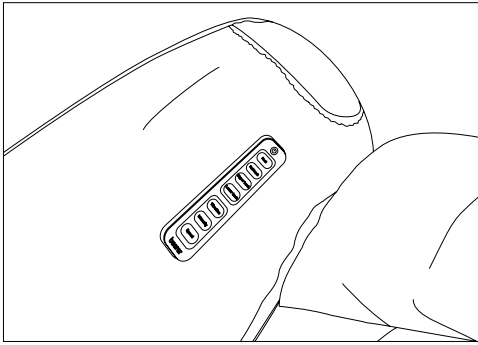- > Connect the power supply unit to the chair using the connector plug
- > Plug the mains plug into a socket
- > Switch on the massage chair using the on/off button  on the control element

The massage chair is now ready for operation.

6. Operation

Warning!

- Before using the device for the first time, read the 'Safety notes' section.
 - ▷ Safety notes, pages 2 and 3


Caution!

- If you experience irritation of any kind, stop treatment immediately.

Before switching on

- Wear comfortable clothes during the massage.
- Remove all objects from your trouser pockets (e.g. keys, pocket knife, mobile phone).
- Assume a comfortable, relaxed position.
- Always sit in the middle of the seat so that your spine is positioned precisely between the massage heads.
- The back cushion can be folded back during the massage.
- The massage should feel pleasant and relaxing at all times. If the massage is painful or uncomfortable, stop the massage or change your sitting position.
- Gently press your back against the back rest.
Now check whether the massage is comfortable.
- Press your weight harder against the back of the chair if this still feels comfortable.
- If the device has not been used for a long period of time, check that the device is functioning correctly before use.

Flat position/upright position

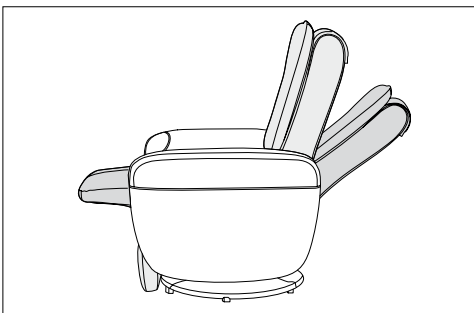The upright position can be adjusted even when the massage function is switched off.


Putting the massage chair in the flat position

- > Sit on the massage chair
- > Pull the straps to release the back rest
- > Push the back rest back with your back
This lowers the back rest and raises the leg rest at the same time.
- > Release the straps as soon as the desired position has been set
This locks the back rest again.

Putting the massage chair in the upright position

- > Pull the straps to release the back rest
- > Push on the leg rest with the legs and release the back rest
This raises the back rest and lowers the leg rest at the same time.
- > Release the straps as soon as the desired position has been set
This locks the back rest again.

Starting the massage

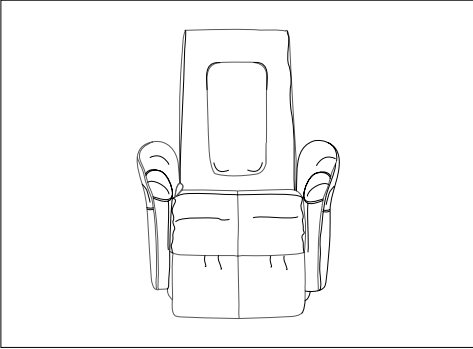- > Plug the mains plug into the socket
- > Fold the back cushion backwards over the back rest
- > Sit on the massage chair
- > Assume the desired upright or flat position
- > Press the on/off button 

The LED button lights up.

- > Press the desired massage button

The massage starts up. The LED button lights up.

The selected programme runs for 15 minutes and then switches off.

The massage chair then switches to standby mode.

Interrupting the programme

- > Press the on/off button 

The massage mechanism then returns to the start position and switches off.

Switching off the massage chair

Once the massage programme is complete, the massage function switches off automatically. Following selection of a new massage area, the massage runs for another 15 minutes.

During the massage, you can switch off the device at any time using the on/off button .

The operating display flashes.

Caution!

Do not use the device for longer than 15 minutes at a time otherwise the device could overheat. Allow the device to cool for at least 30 minutes before using it again.

7. Maintenance and storage

Warning!

- Switch off the device after every use and before cleaning.
- Pull the mains plug out of the socket.


- > Wipe the massage chair from time to time with a dry, lint-free or slightly damp cloth.
- > Every two to three months, wipe the chair using a clean cloth with a small amount of care product for artificial leather.
- > Afterwards, dry the chair with a clean, dry cloth.

Artificial leather is sensitive to grease. Only use care products designated for artificial leather to clean the massage chair.

Important!

- Clean the device using only the methods specified.
- Under no circumstances may liquid enter the device.
- Do not use any solvent-based cleaning products.
- Do not use any harsh cleaning products. These can damage the surface.

Folding mechanism

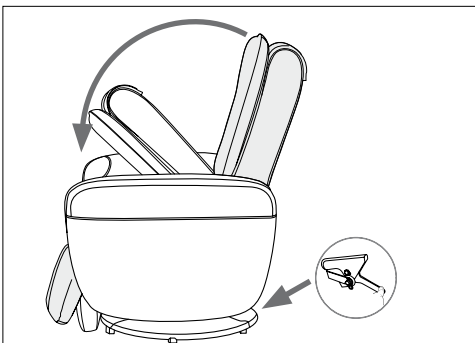The massage chair can be folded away to save space when stored.

- > Pull the mains plug from the socket
- > Release the locking mechanism at the back of the back rest at the bottom
- > Carefully push the back rest forwards on to the seat

If required, lay the back rest cushion back over the back rest before doing so.

8. What if there are problems?

Problem	Cause	Solution
Massage heads are moving slowly.	The load on the massage heads is too great.	> Reduce the contact pressure from your body.
Massage heads are moving in the upper or lower area only.	The massage was only activated for the limited area in question.	> Press the FULL BACK, UPPER BACK or LOWER BACK button to change or extend the area.
Grinding noises during operation.	These sounds are normal and are caused by the operation of the motor or from the massage heads that are rubbing against the cushion.	No action is necessary.
The massage cannot be started.	The mains plug is not connected to the socket.	> Plug the mains plug into the socket.
	The 3D Shiatsu massage was not switched on at the control element.	> Switch on the massage chair using the on/off button  .
The massage function stops when the chair is in operation.	The maximum operation period was reached.	> Switch off the massage chair and allow to cool down for 30 minutes before switching on again.

9. Technical Data

Product description	Massage chair MC3000
Mains voltage/frequency	100-240 V ~ 50-60 Hz
Energy consumption	42 W
Dimensions (in flat position)	Length 140.5 cm, width 76.0 cm, height 99.0 cm
Weight	Net: 30 kg, with packaging 36.5 kg
Maximum load capacity	135 kg

10. Disposal

Observe the local regulations for material disposal. Dispose of the unit in accordance with EC Directive 2002/96/EC – WEEE (Waste Electrical and Electronic Equipment).

If you have any questions, please contact the local authorities responsible for waste disposal.

11. Warranty and service

In the event of a warranty claim or service case, please contact Customer Services in your respective country.

You can find the address in the enclosed address list.

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